

laboratory needs

Flexible

Efficient

Compact

Sample Clot Detection(optional)

Three sensitivity levels of clot detection system to ensure precise sample pipetting.

Automatic Start-up

Analyzer will be always ready for use by running system preparation automatically.

Connectivity

Capable of connecting bidirectionally to LIMS system via ASTM standard.

Automatic Cuvette Check

Cuvettes are thoroughly checked before test to ensure accurate results.

Automatic Rerun

Automatic rerun function enables a true “walk-away” system.

Cuvette

State of art cuvette design makes measurement possible with minimum reaction volume of 150µl.

Reagent Tray

Unique reagent tray system for easy reagent loading.

Wash Unit

Efficient wash system using minimum water and wash solution for eco friendly system.

Clinical Chemistry Analyzer

MODEL

CA-400

Specifications

Clinical Chemistry Analyzer CA-400

Type	Desktop fully automated clinical chemistry analyzer	
Usage	General chemistry as photometric assay Immunology as photometric assay (Latex reagents available)	
Assay type	End point, 2 point end, rate, 2 point rate and electrolytes	
Throughput	400 tests per hour (Up to 560 tests per hour with optional ISE module)	
Incubation time	10 minutes after applying sample for one-reagent assay 5 minutes after applying second reagent for two-reagent assay	
Sample type	Serum, plasma and urine	
Sample input system	Type:	Removable tray with sample tube holder
	Capacity:	Up to 70 samples (handling of STAT available)
	QC sample and Calibrator:	20 positions on designated inner tray
Sampling pipette	Type:	Micro pipette with liquid level detector and clot detector
	Rinsing:	Inside and outside with purified water
Sampling pump	Type:	Micro syringe
	Volume:	2.0 to 35.0 µl (increment by 0.1 µl)
Reagent system	Type:	Removable tray with reagent bottle holder on a turntable
	Capacity:	60 reagent positions on a tray (30 each for 70 ml and 20 ml bottles)
	Inventory:	Calculation of remaining reagent volume available
	Cooling:	Cooling with peltier element (8 to 15 degrees Celsius)
Reagent pipette	Numbers of pipettors:	2
	Type:	Micro pipette with liquid level detector
	Rinsing:	Inside and outside with purified water
	Other function:	Detection of remaining reagent volume available
Reagent pump	Type:	Micro syringe
	R1 Volume:	20 to 350 µl (increment by 1 µl)
	R2 Volume:	20 to 350 µl (increment by 1 µl)
Cuvette system	Number of cuvettes:	90 cuvettes on reaction line
	Type:	Semi-disposable cuvette with onboard washing system
	Material:	Glass
	Volume for measurement: Maximum volume:	Minimum 150 µl 450 µl in total
Reaction system	Type:	Direct heating system
	Temperature:	37.0 degrees Celsius
Detector	Method:	Direct measurement of absorbance in cuvette Bi-chromatic or Mono-chromatic
	Wavelength:	340 nm, 380 nm, 415 nm, 450 nm, 510 nm, 546 nm, 570 nm, 600 nm, 660 nm, 700 nm, 750 nm, 800 nm (12 wavelengths)
	Light source:	Halogen tungsten lamp
Stirring system	Type:	Stick type rotating stirrer
Dimensions	970 mm (W) x 690 mm (D) x 580 mm (H) (Analyzer only)	
Option	ISE module, external tank sensor, clot detection sensor	

Notes Specifications subject to change without notice.

"All brand and product names are registered trademarks,trademarks or service marks of their respective holders."

Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Копомна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37
Пермь (342)205-81-47

Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Саранск (8342)22-96-24
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97
Тверь (4822)63-31-35

Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

https://furuno.nt-rt.ru || fon@nt-rt.ru

Clinical Chemistry Analyzer

MODEL

CA-400

The answer to your laboratory needs

Comprehensive solution to high complex demand

① Highly-Stable Refrigerated Reagent Tray

- ▶ Removable carousel for easy loading, 60 cooled positions to enhance reagent stability.

② Automatic Washing System


- ▶ Effective cleaning method minimizes carry over and contamination.

③ Low Reaction Volume Cuvettes

- ▶ Efficient use of reagent with minimum reaction volume of 150µl.

⑦ Integrated ISE Module (optional)

- ▶ Measurement of Na, K and Cl.
- ▶ Throughput of 240 tests/hour.

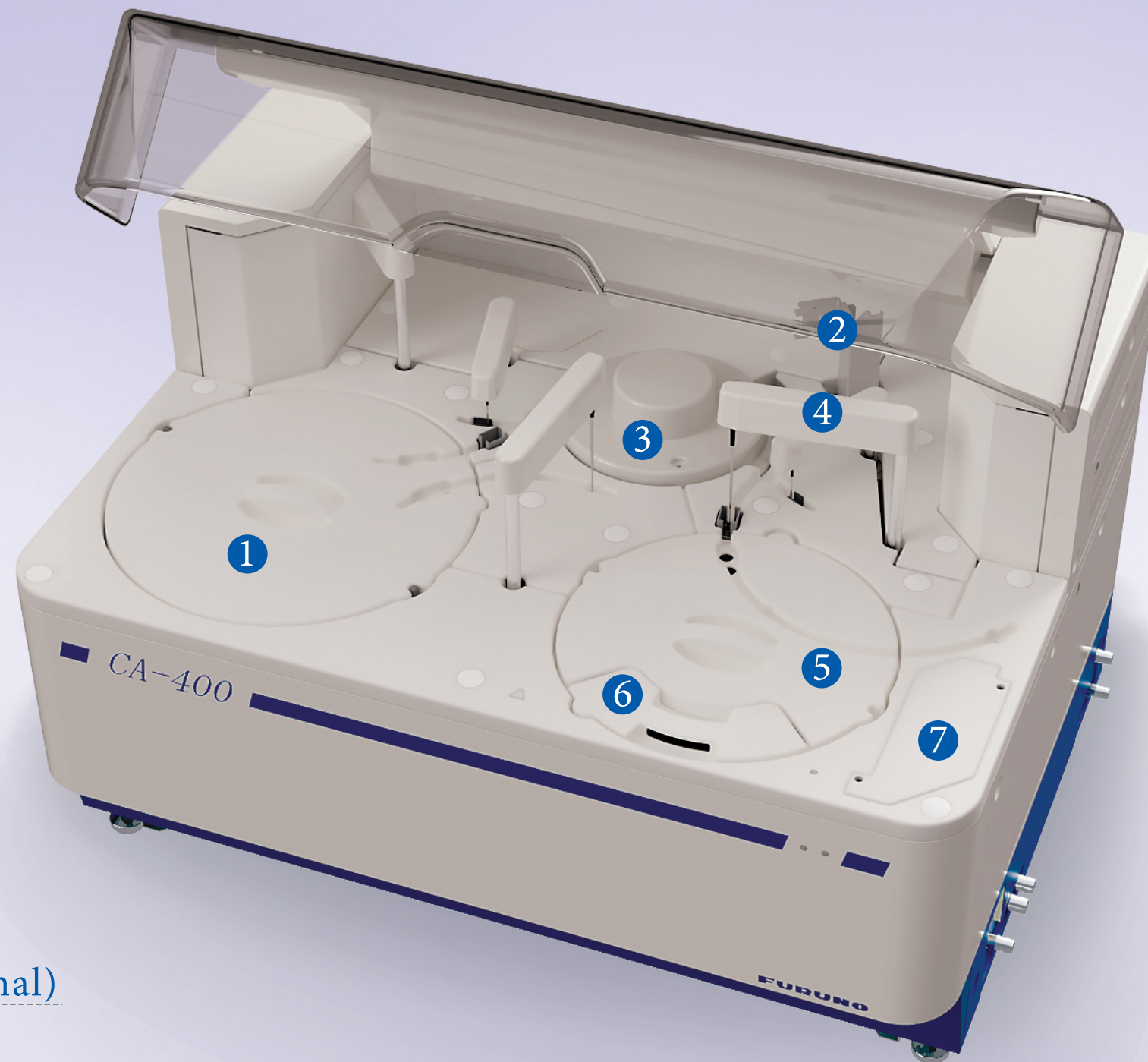

④ High Precision Sample Pipetting System

- ▶ Low sample volume of just 2.0 µl (minimum).
- ▶ Capable of processing various sample types.
- ▶ Liquid level and clash detection.
- ▶ High precision adjustable sample clot detection (optional).

⑤ Multi-Functional Sample Tray

- ▶ Capable of holding primary tubes, sample cups and pediatric cups.


⑥ STAT Sampling

- ▶ Dedicated loading port enables quick process of emergency samples.

User Friendly Interface


- ▶ Intuitive GUI provides fast and simple operation for any level of operator.
- ▶ Key features
 - High operability with easy-to-understand layout
 - Effective use of informative graphics and charts
 - Comprehensive assessment of analyzer condition

Run-Monitor

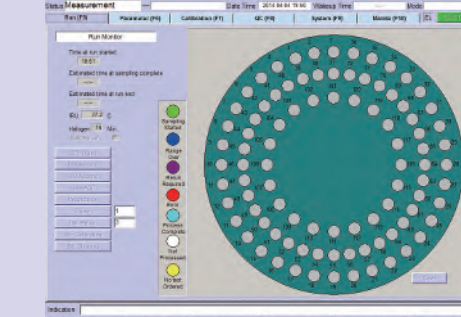


Visually shows sample status.

Quality Control(QC)

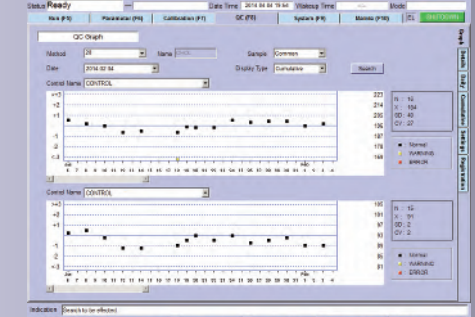


Displays QC results in various formats.

Calibration

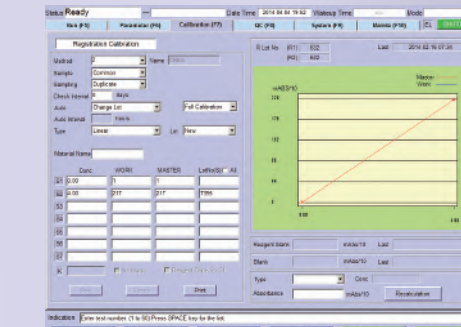


Manages calibration settings and calibration curves on one screen.

Results

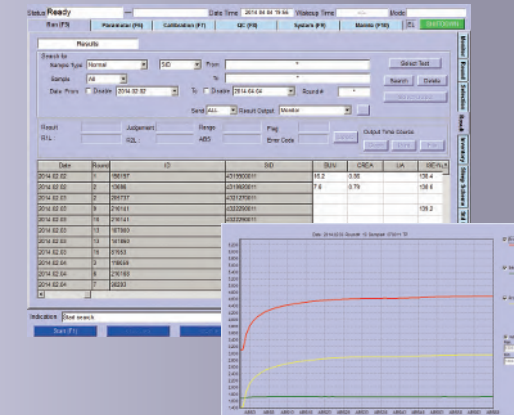

Results can be displayed, printed out, or saved as electronic files.

Dimensions

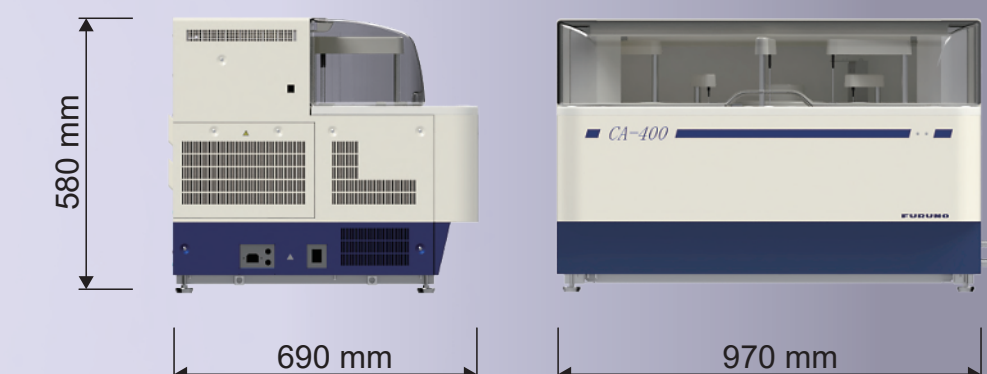

Алматы (7273)495-231	Иваново (4932)77-34-06	Магнитогорск (3519)55-03-13	Ростов-на-Дону (863)308-18-15	Тольятти (8482)63-91-07
Ангарск (3955)60-70-56	Ижевск (3412)26-03-58	Москва (495)268-04-70	Рязань (4912)46-61-64	Томск (3822)98-41-53
Архангельск (8182)63-90-72	Иркутск (395)279-98-46	Мурманск (8152)59-64-93	Самара (846)266-03-16	Тула (4872)33-79-97
Астрахань (8512)99-46-04	Казань (843)206-01-48	Набережные Челны (8552)20-53-41	Санкт-Петербург (812)309-46-40	Тюмень (3452)66-21-18
Барнаул (3852)73-04-60	Калининград (4012)72-03-81	Нижний Новгород (831)429-08-12	Саратов (845)249-38-78	Ульяновск (8422)24-23-59
Белгород (4722)40-23-64	Калуга (4842)92-23-67	Новокузнецк (3843)20-46-81	Севастополь (8692)22-31-93	Улан-Удэ (3012)59-97-51
Благовещенск (4162)22-76-07	Кемерово (3842)65-04-62	Ноябрьск (3496)41-32-12	Саранск (8342)22-96-24	Уфа (347)229-48-12
Брянск (4832)59-03-52	Киров (8332)68-02-04	Новосибирск (383)227-86-73	Симферополь (3652)67-13-56	Хабаровск (4212)92-98-04
Владивосток (423)249-28-31	Коломна (4966)23-41-49	Омск (3812)21-46-40	Смоленск (4812)29-41-54	Челябинск (8352)28-53-07
Владикавказ (8672)28-90-48	Кострома (4942)77-07-48	Орел (4862)44-53-42	Сочи (862)225-72-31	Челябинск (351)202-03-61
Владимир (4922)49-43-18	Краснодар (861)203-40-90	Оренбург (3532)37-68-04	Ставрополь (8652)20-65-13	Череповец (8202)49-02-64
Волгоград (844)278-03-48	Красноярск (391)204-63-61	Пенза (8412)22-31-16	Сургут (3462)77-98-35	Чита (3022)38-34-83
Вологда (8172)26-41-59	Курск (4712)77-13-04	Петрозаводск (8142)55-98-37	Сыктывкар (8212)25-95-17	Якутск (4112)23-90-97
Воронеж (473)204-51-73	Курган (3522)50-90-47	Псков (8112)59-10-37	Тамбов (4752)50-40-97	Ярославль (4852)69-52-93
Екатеринбург (343)384-55-89	Липецк (4742)52-20-81	Пермь (342)205-81-47	Тверь (4822)63-31-35	
Россия +7(495)268-04-70		Казахстан +7(7172)727-132	Киргизия +996(312)96-26-47	

https://furuno.nt-rt.ru || fon@nt-rt.ru